



Wristore™ Distal Radius Fracture Fixator



Visualize, stabilize, mobilize

Wristore™* Distal Radius Fracture Fixator

Today's active adults and young people expect to return to their lifestyles even after surgery. The *Wristore* Distal Radius Fracture Fixator offers an adjustable and stable approach for fractures of the distal radius. By stabilizing fracture fragments near the joint, the *Wristore* Fixator can improve healing and be used to quickly help restore joint function.

A complete external fixation solution

The *Wristore* Fixator is an external fixation approach that:

Provides Extra Stability to Fracture

The *Wristore* Fixator provides extra wrist fixation and increased stability by allowing you to incorporate K-wire fixation to the fixator body. K-wire fixation allows the device to be used in patients with complex wrist fractures.



Is Both Lightweight and Radiolucent

The *Wristore* Fixator is made from a lightweight, high-strength polymer. Its radiolucent housing aids in fracture visualization on x-ray.

* Trademark of Millennium Medical Technologies, Inc.

Multi-Planar Adjustability

The fixator can be easily adjusted using the screws labeled below.

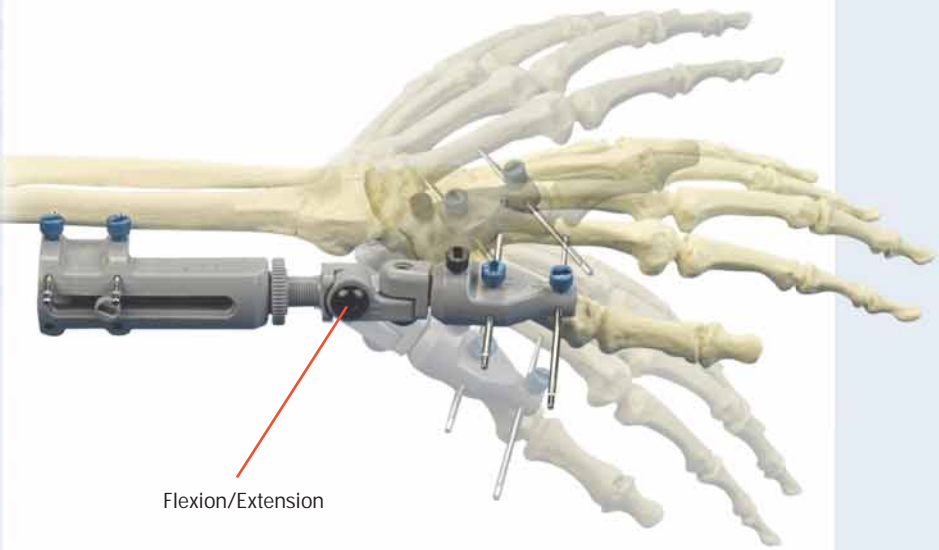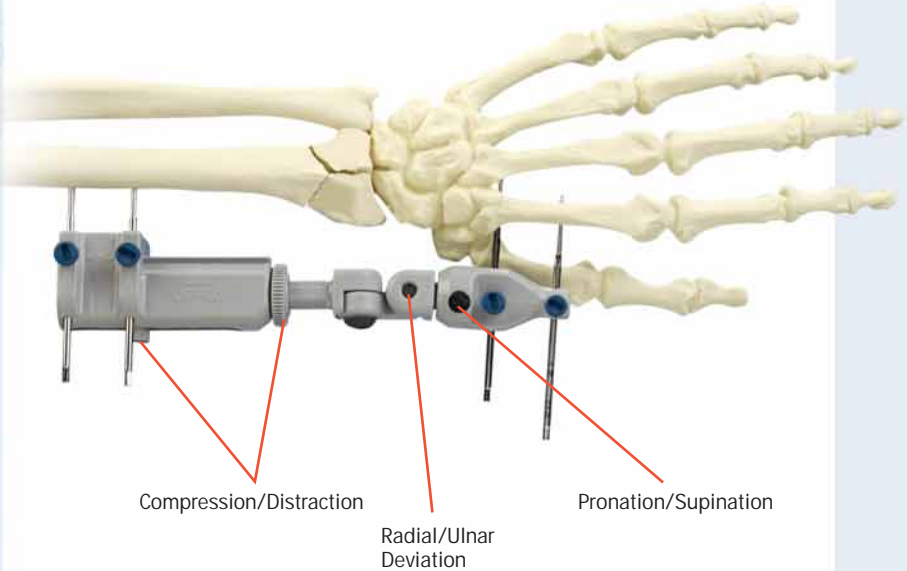



Spanning

Independent Length Adjustment for Controlled Reduction and Stabilization

The *Wristore* Fixator's telescoping body design offers surgeons optimum working length and linear adjustment of the fracture site. This spanning feature is designed to improve fracture stability and joint alignment.

Accurate reduction of the fracture site and direct fracture management and rehabilitation can be achieved through multi-planar adjustability—flexion/extension, pronation/supination, radial/ulnar deviation and compression/distraction.

Complete Sterile Kit

The *Wristore* Fixator surgical kit provides all implants and disposable instruments in one convenient pack, eliminating the need for a sterilization case and extra packaging. A standard OR drill is required for pin and K-wire insertion. To adjust and remove the construct postoperatively, the *Wristore* Pin Inserter/Remover and Wrench/Driver are also available separately.

Surgical Pack

A – Fixator Body

B – Pin Inserter / Remover

C – Five Pin Caps

D – Pin Guide

E – Two 3.0mm Radial Pins

F – Two 2.5mm Metacarpal Pins

G – Dorsal Outrigger

H – Dovetail Slide

I – Five K-wire Caps


J – Four 1.6mm K-wires

K – Two Universal Wire Arms

L – Wrench/Driver

Ordering Information

Prod. No.	Description
00-4600-001-00	<i>Wristore</i> Sterile Pack
00-4600-050-03	<i>Wristore</i> Wrench/Driver
00-4600-052-00	<i>Wristore</i> Pin Inserter/Remover

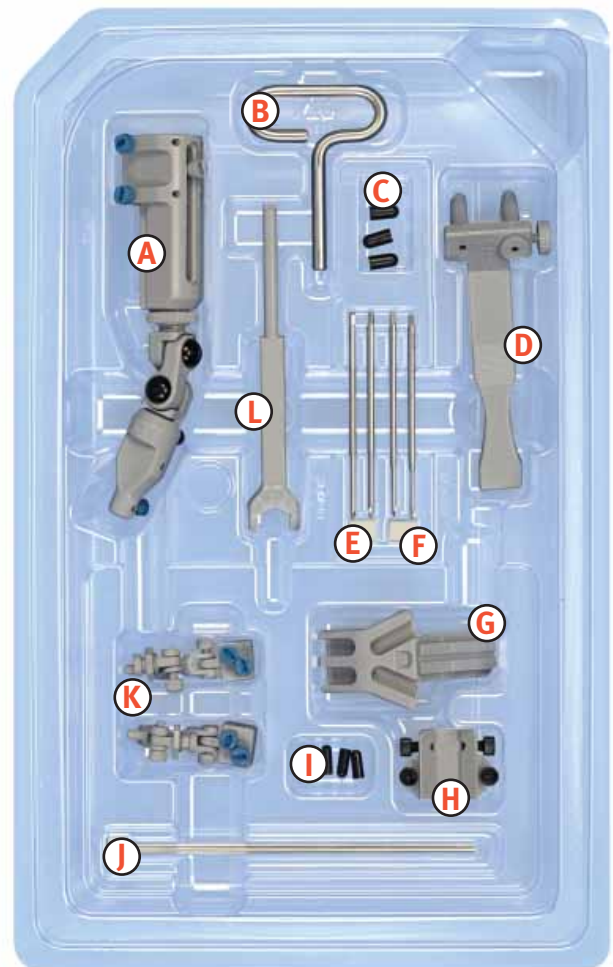

Contact your Zimmer representative or visit us at www.zimmer.com